

# You should Play it safe!

safe!

The RCN Multi-Purpose Shower and  
Commode Chair SCC 250 PPG SB will not let you down ...



safe!

Frame made of high quality  
RCN plastic components

safe!

Swing-out armrest on  
the right ASA

safe!

Supports up to 150 kg

safe!

Safety bar SB

safe!

4 twin-wheel swivel  
casters - all locking

safe!

Anti-tip, slide-out footrest FT

... you should not settle for anything less!

[www.rcn-medizin.com](http://www.rcn-medizin.com)



# Features and Options, to play it



RCN lets you choose the kind of options to make the chair that fits your needs.



Built with high quality RCN plastic components offering longevity and sturdiness. No rust, no sharp edges and high weight capacity provide the ultimate safety for patient and caregiver.




Ø75 mm twin-wheel swivel casters - all locking - keep the chair safe in its place.



The anti-tip, slide-out footrest moves easily back and forth and prevents the chair from tipping forward, even when fully stepped on. The footrest always remains with the chair.



As an additional safety feature the chair can be equipped with an anti-tip frame extension.

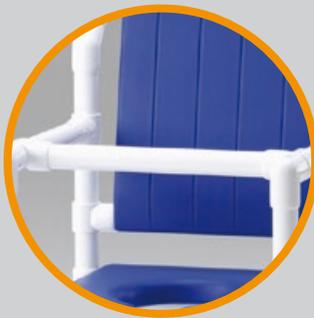

The safety-bar prevents the patient from tipping forward. It is attached to one armrest and simply rotates and clips on to the opposite armrest.



Soft-seat and soft back-rest are made of closed cell PUR foam, both nest on the frame and can be quickly detached for easy cleaning. Front hygiene opening allows better access, a pail comes standard.



The swing-out armrest allows for an easy transfer of the patient on and off the chair. An easy to handle clip secures the armrest.



The adjustable safety-belt is easily attached to the chair frame and keeps the patient from sliding forward.

All components of an RCN chair are guaranteed NO Rust



RCN Medizin- und Rehattechnik GmbH

Hauptstraße 4 · D-55471 Sargenroth  
Tel. +49 6761-9197-0 · Fax +49 6761-9197-70  
info@rcn-medizin.de  
[www.rcn-medizin.com](http://www.rcn-medizin.com)

